



 **Convaid**

Convaid EZ Rider

Lightweight, Compact-Folding, Convenient Fixed-Tilt Buggy

Convaid EZ Rider

The Convaid EZ Rider is a popular upright pediatric buggy that is convenient, lightweight, and compact folding. The EZ Rider is made to order to meet each individual needs. It comes with Convaid's proprietary self-tensioning seating system, Self-Tensio®, which provides increased seating tolerance and can help promote neutral pelvic positioning. It provides growth capability through its seat depth extension system.



Lightweight, Compact Folding, and Convenient



Tension adjustable back upholstery

Adjustable straps allows individual adjustment.



Upholstery with built-in back zippers

Allow for foam inserts for additional padding and support.



Self-Tensio®

Self-tensio is Convaid's unique self-tensioning seating system which provides pelvic positioning similar to a solid seating system in an umbrella folding rehab buggy.



Upholstery with built-in seat zippers

Allow for foam inserts for additional padding and support.

Standard

- 10° Fixed Tilt
- One Piece Height Adjustable Push Handle (EZ 12/14/16)
- Two Piece Push Handles (EZ 18)
- Tension Adjustable Back
- Tool-less Upholstery Removal
- Tool-less Adjustable Seat Depth Growth of 5"
- Zippered Pockets in Seat & Back Upholstery (Cordura Only: EZ 12/14/16)
- 2-Point Positioning Belt (not available in EU)
- Self-Tensio®
- Removable, Height Adjustable Swing-Away Footplates
- 7,5 x 2" Front, 12,5 x 2" Rear Solid Tires (EZ 12/14/16)
- 8" Front, 12,5" Rear Solid Tires (EZ 18)
- Quick Release Rolko Rear Wheels (EZ 12/14/16)
- Foot Operated Wheel Locks (EZ 12/14/16) (Toggle Wheel Locks (EZ 18))
- Silver Powder Coated Frame
- Made in the USA

Options

Sizes

- 12
- 14
- 16
- 18

Seat Depth Tubes

- 51mm seat depth tube extensions
- 127mm seat depth tube extensions (accessory)
- Custom seat depth extension 76mm

Upholstery

- Cordura
- Heavy Duty Cordura (not available in EU)

Tires

- 7,5 x 2" Front, 12,5" Rear Air Knobby Tires
- 12,5" Rear Solid Knobby Tires
- 8 x 2" Front Ribbed, 12,5" Rear Solid or Air Knobby Tires
- 8 x 2" Front Solid Tires

Cordura Upholstery Colours




Princess Pink



Cherry Red



Apple Green

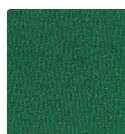

Forest Green



Aqua Blue




Royal Blue



Sassy Purple



Chocolate Brown

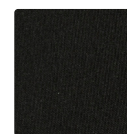

Panther Black

(limited range of colour options in specific countries)

Frame Options



Size 14



Size 18

Accessories



**Extended Headrest Cover
w/ and without windows
(Canopy)**



Headrest Cover (Canopy)



Rain Cover
Canopy Required



Mosquito Net
Canopy Required



**Headrest Extension
Rectangular - Textiene**
(Not available in EU)



**Headrest Extension
Curved - Cordura**



Occi-Headwings



Padded Headwings



Seat Reducer Insert



Full Torso Support Vest



**Soft Adjustable Lateral
Support Single Flap w/
Scoli Strap**



**Soft Adjustable Lateral
Support Double Flap w/
Scoli Strap**

Accessories



H-Harness with padded covers
3-Point Positioning Belt Required



General Use Cushion*¹
(Not available in EU)



Align Cushion - Lateral Pelvis/Lateral Thigh Support*¹
(Not available in EU)



Position Cushion - Medial Thigh Support/Anti Thrust*¹
(Not available in EU)



3-Point Positioning Belt



Lateral Thigh Support (Adductor)



Medial Thigh Support (Abductor)



Height Adjustable Flip-up Armrests



Upper Extremity Support Surface (Clear Tray)



Upper Extremity Support Surface (Foam Tray)[^]



Heavy Duty Reinforced Upholstery - Textiene
(Not available in EU)



Heavy Duty Reinforced Upholstery - Cordura

[^] Crash tested Foam Tray with vinyl cover. Flip-Up Armrests required

*Manufacturer Installation Required

¹Incontinence Cushion Liner available

Accessories



Utility Bag



Travel Bag




Travel Bag for Convertible
Wheels



Calf Panel



Individual Angle
Adjustable Footplates

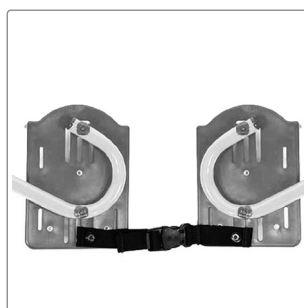

Footplate Securement
Strap



Rear Anti-Tip Tubes (Pair)




Castor Locks
(Standard in EU)



5-Point Harness with Padded
Covers (EZ 12 only)*

Crash Tested, WC19 Compliant
(Not available in EU)

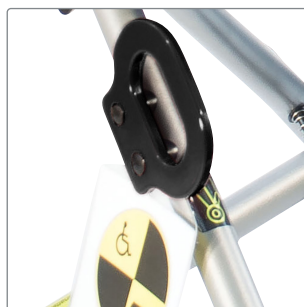

Transit Anchors

*Manufacturer Installation Required

Measurements

EZ Rider measurements		EZ 12	EZ 14	EZ 16	EZ 18
	Seat width (mm)*	305	356	406	457
	Seat depth (mm)*	229-356	279-406	305-432	368-495
	Seat back height (mm)	533	610	660	686
	Seat to floor (mm)	483	508	508	546
	Seat to back angle	100°	100°	100°	100°
	Seat angle	10°	10°	10°	10°
	Overall height (mm)	965	978	1016	940
	Overall length (mm)	914	953	1105	1219
	Overall width (mm)	546	584	660	660
	Folded length (mm)	1130	1168	1219	1156
	Folded width (mm)	356	343	356	368
	Folded height (mm)	445	457	432	394
	Headrest extension (mm)	203	203	203	254
	Shoulder strap height (mm)	279-432	356-508	381-533	406-559
	Seat to footplate (mm)	76-508	76-533	76-559	76-559
	Weight capacity regular/transit (kg)	34/30	45.5/45.5	77/77	114/77
	Weight of chair without front rigging (kg)	11	14	14	13
	Push handle adjustments°	749-1182	737-1194	711-1143	N/A

* Plus or minus 13mm due to fabric
° Adjustable push handle is not available on CX 18 model

Convald-EZRiderBR1201-EUR-REV.00 2020-12-01 | Reserving our rights for printer's error discontinued products. Please find an updated edition on our website © 2020 Etac All rights reserved

For the latest news and contiguously updated product information visit: www.etac.com