



 **Convaid**

Convaid Coaster

Compact-folding, durable, 10°, 20° or 30° adjustable fixed tilt mobility buggy

Coaster

The Coaster is the latest addition to our fixed tilt family. This mobility buggy levels-up the adjustable fixed tilt market by accomplishing simple adjustment between 10°, 20°, or 30° settings, in a lightweight, compactly folded frame that demonstrates the same longevity seen in all Convaid products. The Coaster also offers an adjustable fixed recline across four different positions, providing the optimal comfort and support for all children. Convaid's time-tested accessories from the Cruiser and EZRider are still offered as a part of the Coaster while also presenting new, innovative positioning options.

Standard Features

One piece height adjustable push handle



Tension adjustable back straps



10°, 20° or 30° Adjustable fixed tilt settings

Zippered pockets in seat & back

2-Point positioning belt

Saddle bags

80° Front rigging

Foot operated wheel lock

292mm x 63.5mm Rear solid quick release tyres

Silver powder coated frame

190.5mm x 50.8mm Front solid tyres

Height adjustable swing-away footplates



Made in the USA

Compact-folding, durable, adjustable fixed tilt and recline



Adjustable Fixed Tilt and Recline

With the simple replacement of two bolts, you can adjust the fixed tilt setting to 10°, 20°, or 30° to accommodate the user's needs. For extra positioning, the Coaster offers four fixed recline settings, dependent on the fixed tilt selected.



Zippered Pockets in Seat and Back

The upholstery built in seat and back zippers allows for foam inserts for additional padding and support for the user.



Height and Tension Adjustable Back Positioner

Provides increased lumbar support to encourage proper body positioning and alignment.



Compact Clamshell Folding

Like many commercial strollers, the Coaster can compactly fold into an easily transportable package. The rear wheels are also quick release to offer an even smaller option when loading this buggy.



Continued Convaid Durability

Convaid's reputation-building robustness is once again found in the Coaster, exceeding all industry standards for durability testing.


Add a splash of colour



Sassy purple



Cherry red

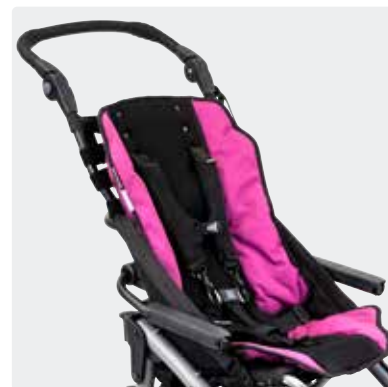


Princess pink



Royal blue



Forest green

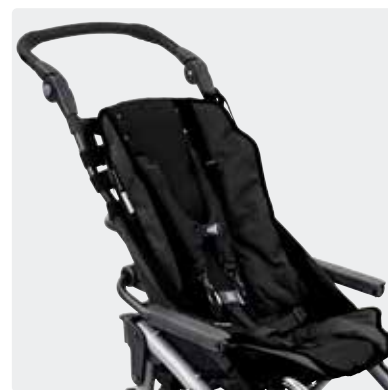

Panther black



Scan to see all product offering



Coaster Scout

The Coaster Scout all-terrain mobility buggy provides families the freedom to truly go anywhere. It comes standard with adjustable fixed seat angle and back angle providing optimal positioning comfort for individuals with mild to moderate trunk positioning needs. The Coaster Scout is designed to withstand rough terrain outdoors with its durable frame and heavy-duty knobby wheels to provide a static grip on more slippery or rocky surfaces. The Coaster Scout also comes with hand-operated riding brakes and the new tension adjustable back positioner for increased safety and support for all adventures.

Standard Features



Made in the USA

Accessories



Mosquito net
Canopy required



Rain cover
Canopy required



**Medical necessity
storage basket***



Utility bag



Headrest cover (canopy)



**Extended Headrest cover
(canopy)**



Headrest extension



Occi-headwings



Padded headwings



**H-harness with padded
covers**



Full torso support vest



**Foam kit for zippered
pockets**



**Height and tension
adjustable back positioner**



**Height adjustable
support tray**



Side panels



3-point positioning belt

*Regular storage basket available



**127mm seat depth
extension tubes**



Lateral thigh support
(adductor)



Medial thigh support
(abductor)



**Angle adjustable
elevating leg rest (90° - 180°)**
(Padded cover available)



Calf panel



**Individual angle adjustable
footplates**



**One piece angle adjustable
flip up footplate**



Full padded footbox
(One-piece footplate required)



Foot positioners



**Height adjustable flip-up
armrests**



Curb tipper



Rear anti-tip tubes



**Caster lock for 50.8mm
wide front wheels**



Transit anchors



Q'straint transit lap belt




















I.V. Pole

New

New

Specifications

Coaster measurements		CTR 12	CTR 14	CTR 16	CTR 18
	Seat width	305mm	356mm	406mm	457mm
	Seat depth	254–381 mm	279–406mm	305–432mm	368–495mm
	Seat back height (10° / 20° / 30°)	559/584/610mm	635/660/686mm	660/686/711mm	711/711/737mm
	Seat to footplate	152–406mm	152–406mm	152–508mm	152–508mm
	Seat to floor (10° / 20° / 30°)	483/508/533mm	483/508/533mm	533/559/584mm	533/559/584mm
	Overall height	1016mm	1041mm	1067mm	1118mm
	Overall length	762mm	787mm	813mm	838mm
	Overall width	559mm	610mm	660mm	711mm
	Folded height	432mm	432mm	432mm	432mm
	Folded length	813mm	813mm	889mm	889mm
	Folded width	559mm	610mm	660mm	711mm
	Headrest extension height	203mm	203mm	203mm	254mm
	Shoulder strap height (10° / 20° / 30°)	267–406/ 292–432 / 318–457mm	330–483 / 356–508 / 381–533mm	356–508 / 381–533 / 406–559mm	419–610 / 445–635 / 470–660mm
	Push handle adjustments	813–1168mm	864–1194mm	864–1270mm	889–1295mm
	Weight capacity regular/transit	34/30kg	45/45kg	77/77kg	77/77kg
	Weight of buggy without front rigging	13kg	13.5kg	14kg	14.5kg
	Fixed tilt and recline settings				
	10° fixed tilt				100°, 105°, 110°, 115° recline options
	20° fixed tilt				90°, 95°, 100°, 105° recline options
	30° fixed tilt				80°, 85°, 90°, 95° recline options

For the intended purpose of the products and latest information, please visit www.etac.com